

Small Sized Storage Series (Square Racks)



Scope of Application:

Widely used within many laboratories, this small-sized storage series features low LN₂ consumption and dual handle design. Stores between 100 and 1100 vials in square racks and cryo boxes.

Innovative Design

- Ideal for Lab Applications
- High Thermal Efficiency
- Lightweight
- Secure
- 30 Litre to 50 Litre Capacity
- 5 Year Vacuum Warranty
- Superinsulation
- Durable

Qingdao Haier Biomedical Co.,Ltd.

No.280 Feng Yuan Road, High-tech Zone,
Qingdao, 266109, P.R. China
Tel: +86-0532-88935593
E-mail: inquiry@haierbiomedical.com
Website: www.haiermedical.com



Haier Biomedical
International



Haier Biomedical
International



@haiermedicalint



Haier Biomedical
International



Haier Biomedical
International

Small Sized Storage Series (Square Racks)

Key Features

- Heavy duty lockable enclosure offers excellent security
- Strong, lightweight aluminium construction
- Dual handles
- Temperature monitoring available
- Roller bases available
- High thermal efficiency
- 5 year vacuum warranty
- Ultra-low LN₂ evaporation loss

Product Advantages

Low LN₂ Consumption and Maximum Storage Efficiency

Haier Biomedical's Small Sized Storage Series (Square Racks) Liquid Nitrogen Storage System features low LN₂ consumption and small storage capacity for laboratory applications. The lightweight containers have a smaller footprint whilst providing maximum efficiency of cryogenic storage capacity.

Advanced Vacuum Technology and Super-insulation Technology


Advanced vacuum and insulation technologies provide ultra-low liquid nitrogen evaporation loss and ensure up to three months cryogenic storage.

Five Models Available

Small Sized Storage Series (Square Racks) offers Five models with different capacity from 100, 600, 750, 900 to 1100 standard 2 ml cryovials. All models are supplied with stainless steel racks and cryogenic boxes.



Accessories

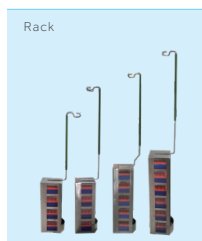

Small Sized Storage Series (Square Racks)

Technical Parameters



Model	YDS-10-125	YDS-30-125	YDS-35-125	YDS-47-127	YDS-50B-125
Maximum Storage Capacity					
No. of Rack	1	6	6	6	6
1.2&2 ml Cryovials (25 per box)	100	600	750	1110	900
No. of Boxes in (Rack)	4	4	5	5	6
Performance					
Volume of LN ₂ (L)	10	31.5	35.5	47	50
Static Evaporation* (L/Day)	0.42	0.35	0.36	0.36	0.45
Static Holding Time** (Day)	24	90	97	130	110
Dimensions					
Opening Diameter (mm)	125	125	125	127	125
Height (mm)	625	704	748	718	818
Outside Diameter (mm)	300	462	462	508	462
Empty Weight (kg)	6.3	13	14.5	18.2	17.3
Full Weight (kg)	15.1	38.9	43.2	60.9	62.4
Shipping Weight (kg)	9.6	19.3	21.9	30.7	30.3

* Static evaporation and static holding time are theoretical values. Actual evaporation and holding time are affected by usage, atmospheric conditions and manufacturing tolerances.

** Static Holding Time: the amount of days it takes for all of the LN to evaporate out after container is filled with liquid nitrogen to pre-cool and reach thermal equilibrium and the liquid nitrogen filling rate reaches 100% under standard atmospheric condition.

Product List



* All models with F-suffix are supplied with stainless steel racks and cryo boxes.

Model	Product	Model	Description and Accessories
YDS-10-125	Liquid Nitrogen Container	YDS-10-125	No canister; no protective jacket; no lockable cap; equipped with braces for carrying
	Liquid Nitrogen Container	YDS-10-125-F	Include 1 square rack with 4 layers and 5*5 cryobox; lockable cap; no protective jacket
	Rack	FTT-25-4	Used to store 2ml cryo boxes (5*5), 4 layers in total, 5*5 cryo boxes included, suitable for YDS-30-125/YDS-10-125 (plastic cryo boxes included)
	Canister (TT97)	TT10H-125	Φ97 single tier canister, no goblet, suitable for YDS-10-125, YDH-10-125
	Cork (125)	GS10-125	Suitable for YDS-10-125, YDH-10-125, YDS-50B-125
	Lockable Cap	SG-125	Suitable for container with 125mm neck opening
YDS-30-125	Liquid Nitrogen Container	YDS-30-125-F	Includes 6 racks with 4 layers, 5*5 cryoboxes, lockable cap and protective jacket, roller base is optional.
	Liquid Nitrogen Container	YDS-30-125	Includes protective jacket only, roller base is optional
	Rack	FTT-25-4	Used to store 2ml cryo boxes (5*5), 4 layers in total, 5*5 cryo boxes included, suitable for YDS-30-125/YDS-10-125 (plastic cryo boxes included)
	Canister (TT97)	TT30-125	Φ97 single tier canister, no goblet, suitable for YDS-30-125
	Roller Base	YSC-S (2019)	Suitable for all 30L / 35L / 50L aluminium tank, equipped with 5 universal wheels
	Smart Cap	LT-125-370	Used to monitor level and temperature for containers with 125mm neck opening, data can be transmitted to cloud platform, suitable for YDS-20-125/YDS-30-125
	Cork (125)	GS20-125	Suitable for YDS-30-125, YDS-35-125, YDS-47-127, YDS-35B-125
	Lockable Cap	SG-125	Suitable for container with 125mm neck opening

Small Sized Storage Series (Square Racks)

Product List



Model	Product	Model	Description and Accessories
YDS-35-125	Liquid Nitrogen Container	YDS-35-125-6YS	Includes 6 canisters (height 276mm), lockable cap and protective jacket, roller base is optional
	Liquid Nitrogen Container	YDS-35-125-F	Includes 6 racks with 5 layers, 5*5 cryoboxes, lockable cap and protective jacket, roller base is optional.
	Liquid Nitrogen Container	YDS-35-125	Includes protective jacket only, roller base is optional
	Rack	FTT-25-5	Used to store 2ml cryo boxes (5*5), 5 layers in total, 5*5 cryo boxes included, suitable for YDS-35-125, YDS-47-127 (plastic cryo boxes included)
	Canister (TT97)	TT35-125	Φ97 single tier canister, no goblet, suitable for YDS-35-125
	Canister (TTS97)	TTS35-125	Φ97 double-tier canister, 2 goblets included, suitable for YDS-35-125
	Roller Base	YSC-S (2019)	Suitable for all 30L / 35L / 50L aluminium tank, equipped with 5 universal wheels
	Smart Cap	LT-125-400	Used to monitor level and temperature for containers with 125mm neck opening, data can be transmitted to cloud platform, suitable for YDS-25-125 / YDS-35-125 / competitors' products
	Cork (125)	GS20-125	Suitable for YDS-30-125, YDS-35-125, YDS-47-127, YDS-35B-125
	Lockable Cap	SG-125	Suitable for container with 125mm neck opening
YDS-47-127	Liquid Nitrogen Container	YDS-47-127-6YS	Includes 6 round rack with 5 layers (37 holes) and lockable cap, no protective jacket, roller base is optional
	Liquid Nitrogen Container	YDS-47-127	Includes 10 canisters (Φ 72 / 276mm) and lockable cap, no protective jacket, roller base is optional
	Round Rack (Y98-37)	TT47-127-5	Φ98 round rack used to store 37 vials (2ml), 5 tiers in total, suitable for YDS-47-127-6YS
	Canister (TTS72)	TTS47-127-10T	Φ72, height 276mm double-tier canister, 2 goblets included, suitable for YDS-47-127
	Roller Base	YSC-47L	Suitable for all 47L Aluminium tank, equipped with 5 universal wheels
	Smart Cap	LT-127-370-10	Used to monitor level and temperature for containers with 125mm neck opening, data can be transmitted to cloud platform, suitable for YDS-47-127 / competitors' products
	Smart Cap	LT-127-400	Used to monitor level and temperature for containers with 125mm neck opening, data can be transmitted to cloud platform, suitable for YDS-47-127-6YS / YDS-47-127 / competitors' products
	Cork (125)	GS20-125	Suitable for YDS-30-125, YDS-35-125, YDS-47-127, YDS-35B-125
	Cork (125)	GS47-127-10T	Suitable for YDS-47-127
	Lockable Cap	SG-125	Suitable for container with 125mm neck opening
YDS-50B-125	Liquid Nitrogen Container	YDS-50B-125-F	Includes 6 racks with 6 layers, 5*5 cryoboxes, lockable cap and protective jacket, roller base is optional.
	Liquid Nitrogen Container	YDS-50B-125-S	Includes 6 canisters (Φ97 / 276mm), lockable cap and protective jacket, roller base is optional
	Liquid Nitrogen Container	YDS-50B-125	Includes protective jacket only, roller base is optional
	Rack	FTT-25-6	Used to store 2ml cryo boxes (5*5), 6 layers in total, 5*5 cryo boxes included, suitable for YDS-50B-125 (plastic cryo boxes included)
	Canister (TT97)	TT50-125	Φ97 single tier canister, no goblet, suitable for YDS-50B-125
	Canister (TTS97)	TTS50-125	Φ97 double-tier canister, 2 goblets included, suitable for YDS-50B-125
	Roller Base	YSC-S (2019)	Suitable for all 30L / 35L / 50L aluminium tank, equipped with 5 universal wheels
	Smart Cap	LT-125-500	Used to monitor level and temperature for containers with 125mm neck opening, data can be transmitted to cloud platform, suitable for YDS-50-125
	Cork (125)	GS10-125	Suitable for YDS-10-125, YDH-10-125, YDS-50B-125
	Lockable Cap	SG-125	Suitable for container with 125mm neck opening